




LG Pra.L Derma LED Mask

Model name
BWJ1V




1903_Rev01

Charge

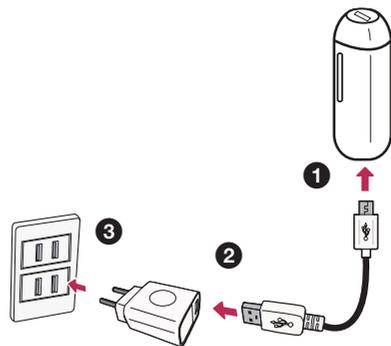

Power on / off

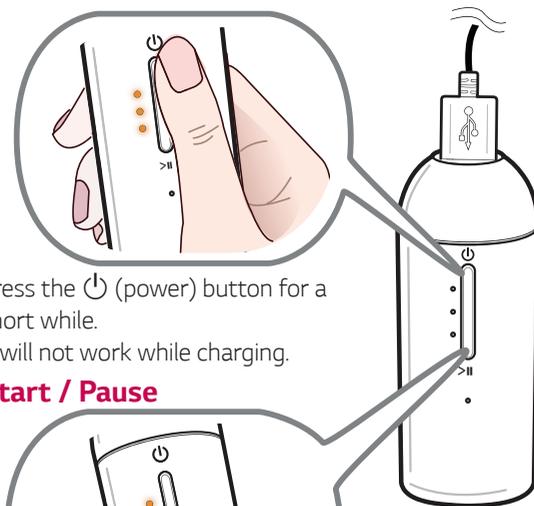


Press the  (power) button for a short while.
It will not work while charging.

Start / Pause



Press the  (Start/Pause) button for a short while.

Set voice guide function

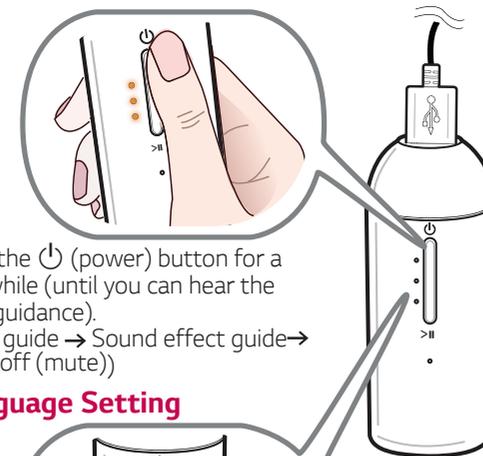

Press the  (power) button for a long while (until you can hear the voice guidance).
(Voice guide → Sound effect guide → Guide off (mute))

Language Setting



Press the  (power) button and the  (Start/Pause) button at the same time for a long while (until you can hear the voice guidance).
(Chinese → English)

Safety cautions

- When using this device, depending on individual skin condition, it may cause a variety of skin irritation. If such irritation persists, please reduce the number of use.
- In the following cases, it is recommended to use the device after consultation with a specialist.
 - Allergies skin, acne, dermatitis, skin treatment, cosmetic surgery or hair loss surgery
 - Wearing a cardiac pacemakers, device implantation or a metal prosthetic device
- As a facial treatments only, it is not recommended to use for a body.
- Do not make direct eye contact with LED light in operation. Excessive exposure to the light can cause eye damage.
- Recommend using the product once a day. Using more than twice a day may cause burning sensation on skin.
- Before using the device, if you are currently undergoing dermatology related treatment or have skin problems, consult your dermatologist and then use it.
- Pregnancy, menstruation, lactating women are not recommended for use. (Unstable hormonal balance may cause temporary skin trouble.)
- Do not use the device for purposes other than those set out in this guideline.
- Please pay attention to using cosmetics that contain the following ingredients as it may reduce the effect. (For example: denatured alcohol, sunscreen, retinol, oil)




How to Use




Wash the face.



Apply the appropriate skin / lotion to your skin.

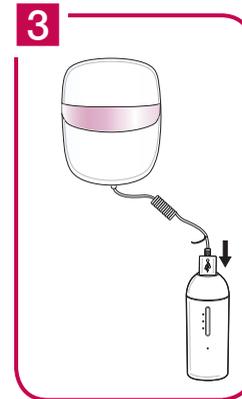

Connect the cable of Derma LED Mask to the port of the controller.

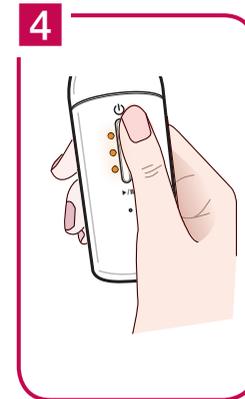

Turn on by pressing the  (Power) button.



Wear the Derma LED Mask as if wearing eye glasses.



Press the  (Start/Pause) button to begin skincare.

① If it does not work after wearing the mask, please let your skin close adhere to the proximity sensor in the middle of the wearing part.