




LG Pra.L Galvanic Ion Booster


Model name
BBJ1 BBJ2



Charge

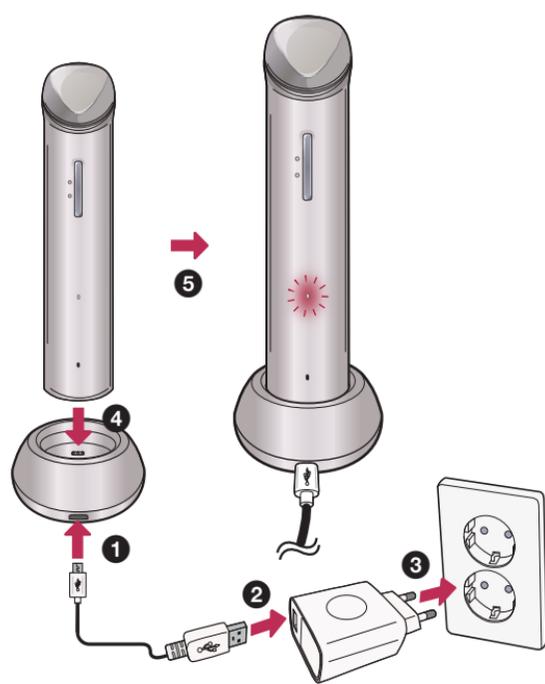


Power on / off

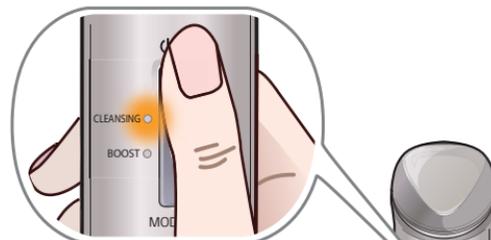


Press the (power) button for a short while.

Mode Selection (CLEANSING / BOOST)

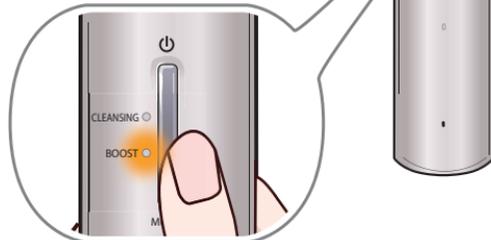

Press the MODE button for a short while.
(CLEANSING → BOOST)

Set voice guide function

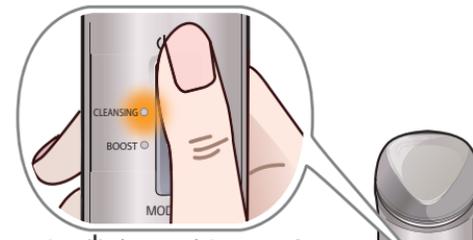

Press the (power) button for a long while (until the guide).
(Voice guide → Sound effect guide → Guide off (mute))

Language Setting



Press the (power) button and the MODE button at the same time for a long while (until the guide).
(Korean → English → Chinese)

Safety cautions

- When using this device, depending on individual skin condition, it may cause a variety of skin irritation. If such irritation persists, please reduce the number of use.
- In the following cases, it is recommended to use the device after consultation with a specialist.
 - Allergies skin, acne, dermatitis, skin treatment, cosmetic surgery or hair loss surgery
 - Wearing a cardiac pacemakers, device implantation or a metal prosthetic device.
- As a facial treatments only, it is not recommended to use for a body.
- Do not expose the LED light of the UV head protection cap directly. Excessive exposure to the light can cause eye and skin damage.
- Recommend using the product once a day by mode. Using more than twice a day may cause burning sensation on skin.
- Before using the device, if you are currently undergoing dermatology related treatment or have skin problems, consult your dermatologist and then use it.
- Pregnancy, menstruation, lactating women are not recommended for use. (Unstable hormonal balance may cause temporary skin trouble.)
- Do not use the device for purposes other than those set out in this guideline.
- Please pay attention to using cosmetics that contain the following ingredients as it may reduce the effect. (For example: denatured alcohol, sunscreen, retinol, oil)

Cleansing Mode (CLEANSING)

※ To make your skin cleaner, it is effective to use this mode before using Absorption Promotion Mode.

Main functions: Heat, Ion



1 Remove makeup and wash the face.




2 Apply moisture gel.



3 Turn on the power.



4 Press the MODE button and select the CLEANING Mode.

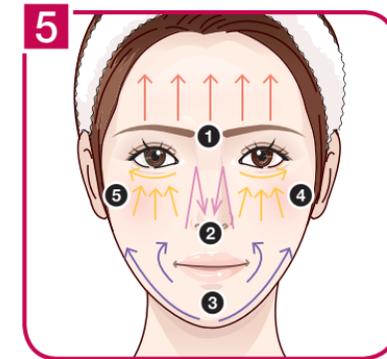

According to the Voice Guide function, cleanse each area of the face in the order in the direction of an arrow.

- ➔ 1 Forehead (45 seconds)
- ➔ 2 Nose (45 seconds)
- ➔ 3 Chin (30 seconds)
- ➔ 4 Left cheek (30 seconds)
- ➔ 5 Right cheek (30 seconds)



6 Rinse your face sufficiently with tepid water after cleansing.

Absorption Promotion Mode (BOOST)

※ The complex function of the mode helps the penetration of active substances contained in cosmetics.

Main functions: Heat, Ion, Ultrasound



1 Remove makeup and wash the face.




2 Apply cosmetics (eg, Vitamin C Serum).



3 Turn on the power.



4 Press the MODE button and select the BOOST (absorption promotion) Mode.

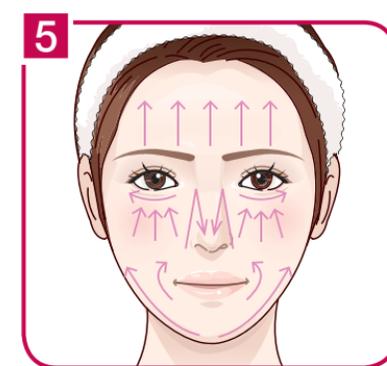

Take care of each desired area for 3 minutes in the direction of an arrow.