



## LG Pra.L Derma LED Mask

Model name  
BWJ1 BWJ2




MBM66277328


1808\_Rev02

### Charge


Charge  
on the  
cradle

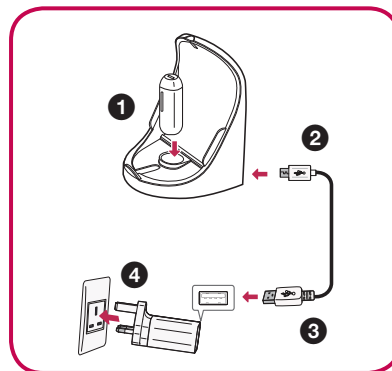

Charge  
directly  
on the  
controller

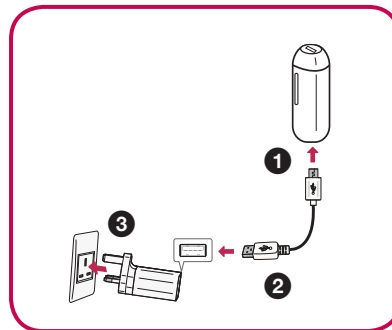


### Power on / off

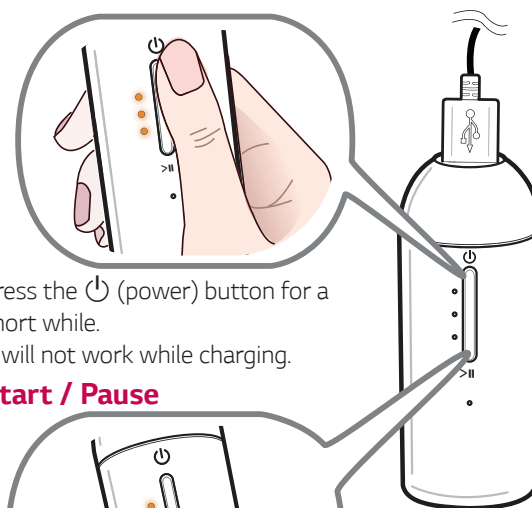


Press the ⏻ (power) button for a short while.  
It will not work while charging.

### Start / Pause

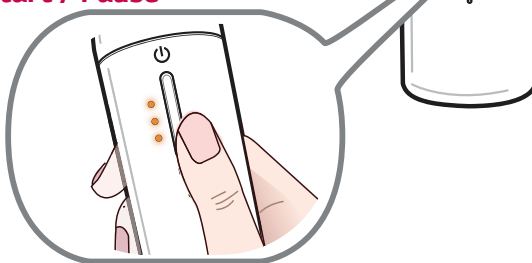

Press the ⏮/⏭ (Start/Pause) button for a short while.

### Set voice guide function

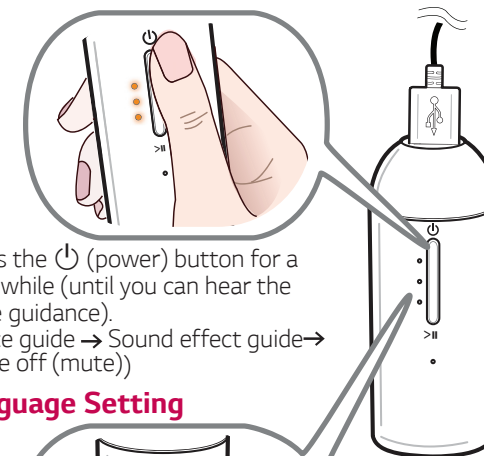

Press the 🔊 (power) button for a long while (until you can hear the voice guidance).  
(Voice guide → Sound effect guide → Guide off (mute))

### Language Setting



Press the 🔊 (power) button and the 🗣 (Start/Pause) button at the same time for a long while (until you can hear the voice guidance).  
(Hong Kong Chinese → English → Chinese)

### Safety cautions

- Stop using the device if there is smoke or an odor. In such cases turn off the power immediately, pull out the charging adapter and contact the service center.
- Always use the provided charging adapter.
- Damaged battery may cause explosion.
- Do not disassemble, repair, or remodel this product.
- Handle the power plug with dry hands.
- Do not make direct eye contact with LED light when devices with LED is in operation. Excessive exposure to the light can cause eye damage.
- Recommend using the product once a day. Using more than twice a day may cause burning sensation on skin.
- Before using the device, if you are currently undergoing dermatology related treatment or have skin problems, consult your dermatologist and then use it.
- Pregnancy, menstruation, lactating women are not recommended for use. (Unstable hormonal balance may cause temporary skin trouble.)
- Do not use the device for purposes other than those set out in this guideline.




## How to Use




Wash the face.




Apply basic skincare products.

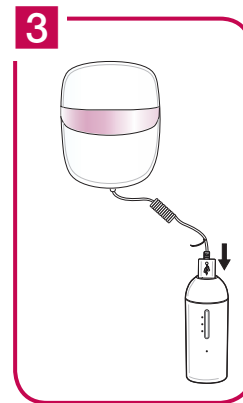

Connect the cable of Derma LED Mask to the port of the controller.

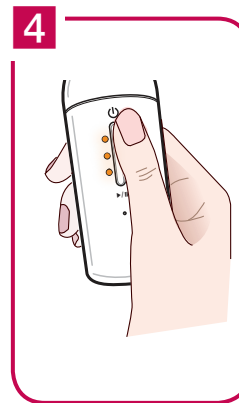

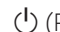
Turn on by pressing the  (Power) button.



Wear the Derma LED Mask as if wearing eye glasses.



Press the  (Start/Pause) button to begin skincare.

① If it does not work after wearing the mask, please let your skin close adhere to the proximity sensor in the middle of the wearing part.