


Instructions for Use

This device is a device used for the treatment of hair loss for androgenic alopecia according to the following categories.

- Norwood Hamilton Classification: Male with androgenic alopecia from IIa to V
- Ludwig Classification: Female with androgenic alopecia from I to II
- BASP specific type Classification: Androgenic alopecia from V1 to V2 and F1 to F2
- Fitzpatrick Classification of Skin Phototypes: Skin type from I to IV


Norwood Hamilton Classification

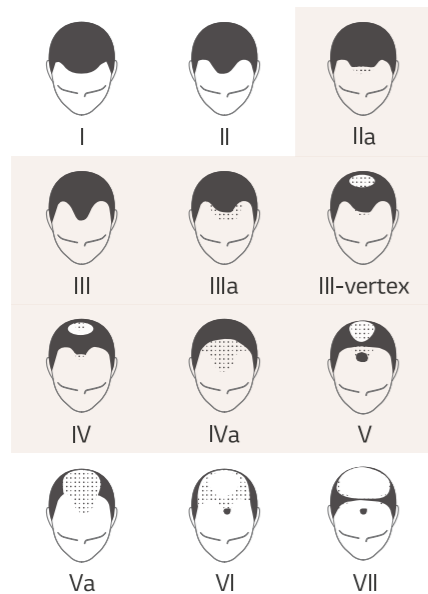

Ludwig Classification



BASP Specific Type Classification

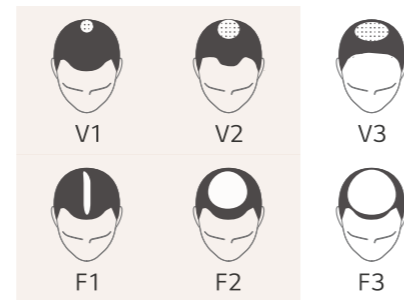

Fitzpatrick Classification of Skin Phototypes



* This device is a medical device.



COV37374102
QSG-01 (Rev. 07)
2022. 12.
2212_Rev03

LG Pra.L MEDI HAIR



For more information, refer to the provided owner's manual or the online owner's manual of LG web site. You can quickly access the user manual by scanning the QR code with a smart device. This device is a medical device.


HGN1

Precautions for use

- Recommended number of use: 3 times a week
 - Do not use this device for 2 days in a row.
 - Using this device for longer than the recommended duration may cause a burning sensation on scalp.
- Do not use this device if you have any injury or disorder on scalp or any photosensitivity reaction such as solar urticaria.
- Device to be used as a scalp treatments only, it is not recommended for use on body.
- Device use not recommended for pregnant, menstruating or lactating women.
- Depending on individual scalp condition, the device may cause a variety of scalp irritations. If such irritation persists, please reduce the number of uses.
- In the following cases, please use the device after consultation with a specialist.
 - If you are in the recovery period after treatment related to scalp disease, forehead and hair
 - If you have any ophthalmologic diseases or recently had an ophthalmologic operation
 - If your skin is sensitive or allergic
 - Those with cancer, acute disease, high fever, cardiovascular disorder, endocrinopathy or those with any physical disorder
 - Those using implanted medical device or prosthetics under/ around scalp exposed to laser/LED light
 - Those suffering from itchiness, eye pain, dizziness, headache, dermalgia, infection, skin disorder, erythema, pimples, and asteatosis
 - Do not use the product for a week if you have received hair treatment such as perm, henna, dyeing, and bleaching
- Do not use the device for purposes other than those set out in this guide.

Getting Ready for Care

- 1 It is recommended to maintain the hair and scalp clean when using the product in order to avoid any hindrance to the laser/LED irradiation. (When using the product, remove any hair-styling product or wig from your head.)
- 2 If you want to use the product after taking a shower or washing your hair, make sure that you dry the hair and scalp sufficiently before using the product. (It is normal for you to feel burning sensation or sweat on the scalp while using the product.)
- 3 Put on the main unit in a comfortable posture.

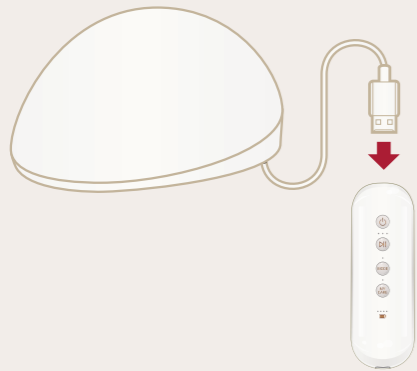


Connect the Main Unit with the Controller

Connect the LG Pra.L MEDI HAIR main unit with the controller.

- Before using the device, make sure the controller is fully charged. Use after the device is fully charged.
- ※ Be sure to use the power adapter included in the package.

How to Use

Recommend using 3 times a week

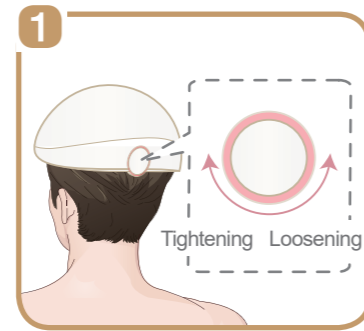


Put on the main unit and adjust it to your head size with the band adjustment dial on the back.

※ Rear cushion pad (aid pad for adjusting the size)



- If the cushion pad doesn't fit your head because it is too large, replace the pad with a thicker pad additionally provided as an accessory.

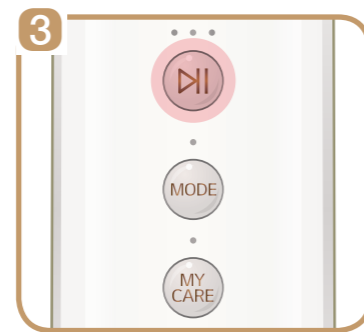

Select the desired care mode and start the care by pressing the Start/Pause button.


- MODE Setting: Press the MODE button for a short time to change the mode. (Total Care → Front/Upper hair care → Upper hair/Top head care)
- MY CARE: Recommended care mode according to user's hair loss condition

※ The mode cannot be changed once the care starts.



Turn on the device by pressing the Power button on the controller.

Using the MY CARE function

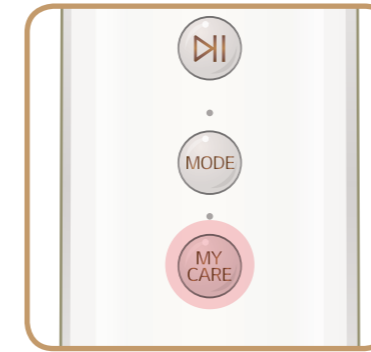

With power on, press the MY CARE button for a short time and press the Start/Pause button for a short time.

The MY CARE function identifies your hair loss type, recommends the care mode for you, and save the information.

※ Press the Start/Pause button to pause the care temporarily. If the care pauses for more than 5 minutes, the power turns off automatically.


※ If you want to re-measure the hair loss condition of your saved MY CARE, press and hold the MY CARE button for 5 seconds.

Care Mode	Status
Total Care (27 minutes)	<p>9 minutes (Front/Upper hair) → 9 minutes (Upper hair/Top head) → 9 minutes (Front hair/Top head)</p> <p>All parts are cared comprehensively. This mode is for those who want total care or who have poor hair thickness and density of the entire hair.</p>
Front/Upper hair care (18 minutes)	<p>Front/Upper hair are cared. This mode is for those who have the front line pushed backwards, who have poor hair thickness and density of the upper hair.</p>
Upper hair/Top head care (18 minutes)	<p>Upper hair/Top head are cared. This mode is for those who have poor hair thickness and density of the upper hair and top head.</p>

※ For more information on the care modes, refer to the "Using the MODE" section in the user manual.


※ For the hair loss condition depending on men and women, refer to 'Instructions for Use' in the back side of this quick setup guide or the owner's manual.

Using the MODE

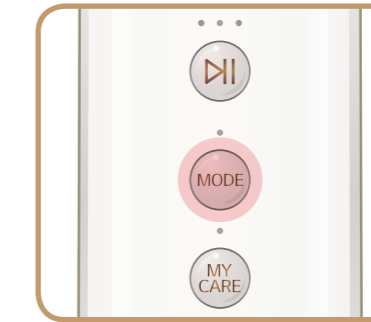

With power on, press the MODE button for a short time to cycle through each mode: (Total Care → Front/Upper hair care → Upper hair/Top head care) Select the desired care mode and start the care by pressing the Start/Pause button for a short time.

Voice Assistance Setting

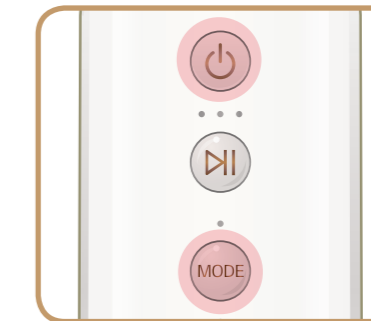

With power on, press and hold both the Power button and the MODE button simultaneously for about 3 seconds to cycle through each mode: (Voice assistance → Sound effect guide → Turn off guide (mute))

Setting up Language



With power on, press and hold both the Power button and the Start/Pause button simultaneously for about 3 seconds to change each language.

※ When turning off the device, press the Power button to turn it off. Then the setting information such as care mode, language, and voice assistance is automatically saved.