



LG Pra.L

SIMPLE MANUAL

Skin Tightening Device

Please read this manual carefully before operating your set and retain it for future reference.

Model
BLP1



Detailed instructions can be found in the WEB manual from LG Electronics official website. You can quickly access the user manual by scanning the QR code with a smart device.



COV37314803
2210_Rev06

www.lg.com

Product Features

This device uses an advanced technology of radiofrequency / microcurrent for elasticity care, and another advanced technology of sonophoresis / iontophoresis for skincare absorption.

※ Radiofrequency: A technology which helps improve elasticity by generating deep heat inside the skin tissue

※ Microcurrent: A technology which improves elasticity by stimulating skin tissues and muscles

- What is microcurrent? Microcurrents are low-level currents under 1000 μ A, similar to bioelectricity flowing in the body

※ Sonophoresis: A technology which delivers skincare ingredients into the skin through ultrasonic waves

※ Iontophoresis : A technology of penetrating the effective components of skincare by applying a potential difference to the skin

Intensive Care

This device helps improve skin elasticity and promote skincare absorption simultaneously as “one product” from “one surface.”


Eye Care

This device has a dedicated Eye Care Head that helps promote younger-looking eyes using hypoallergenic complex elasticity care technology.


Cooling Care

Cooling Care helps to minimize the appearance of pores on the skin.

Components

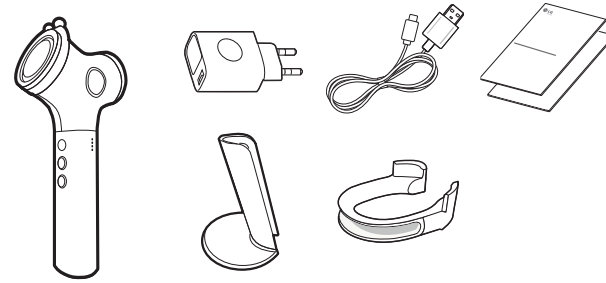


Unit Description

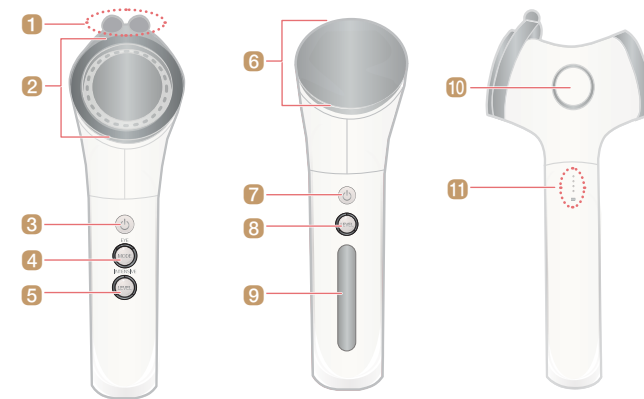

- 1 Eye Care Head
- 2 Intensive Care Head
- 3 Intensive / Eye Care Power (ON) Button
- 4 Intensive / Eye Care MODE (MODE) Button, MODE Status Light
- 5 Intensive / Eye Care LEVEL (LEVEL) Button, LEVEL Status Light
- 6 Cooling Care Head
- 7 Cooling Care Power (ON) Button
- 8 Cooling Care LEVEL (LEVEL) Button, LEVEL Status Light
- 9 Skin Contact Electrode
- 10 Ventilation Hole
- 11 Battery Status Light

Charge

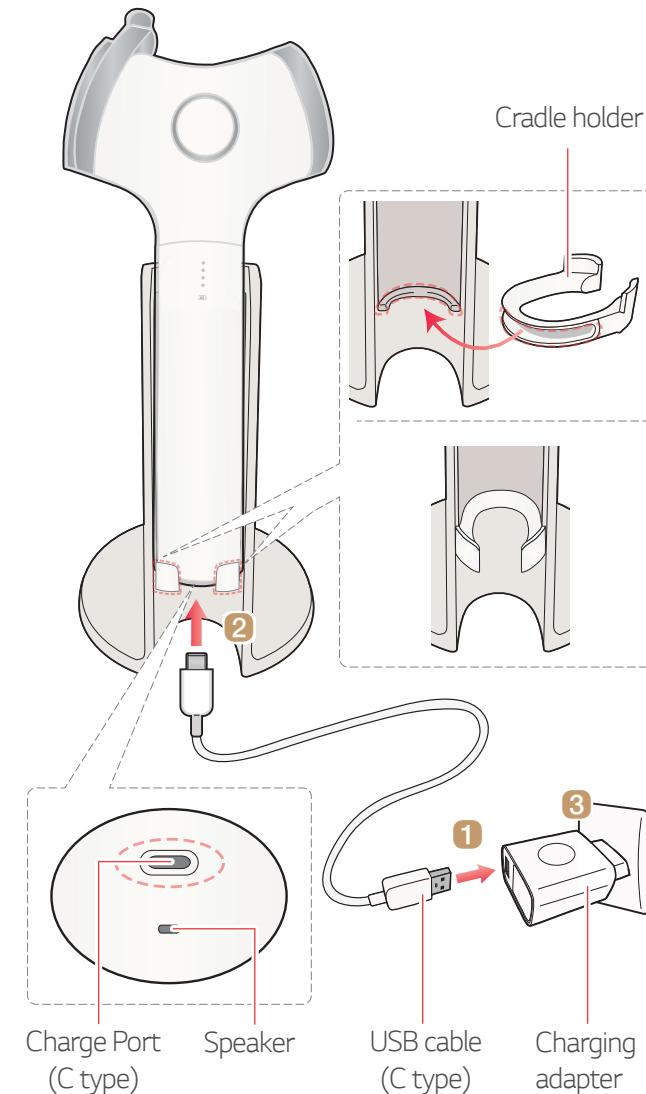

! This unit cannot be used while charging.

Product specification

General	
Power (Power adapter)	Model : ADS-12BA-06Y 05010EPK C Manufacturer : Shenzhen Honor Electronic Co., Ltd. Input : 100 - 240 V ~ 50 / 60 Hz Output : 5 V $\overline{\text{---}}$ 2.0 A
Dimensions (L x H x W)	Main Body : Approximately 42.3 mm x 215 mm x 33.6 mm Cradle : Approximately 115.8 mm x 148.5 mm x 11.3 mm
Operating temperature	5 °C - 35 °C
Operating humidity	0 % - 90 %
Recommended number of use	3 times a week (based on the Intensive Care Mode)
	• After being fully charged, it can be used for about 60 minutes based on Level 3 of Intensive Care. (Internal test result) • The hours of use can vary depending on the battery warranty period, surrounding temperature and conditions of use.
Time to full charge	Maximum of 3 hours
	• Based on 5 V, 2.0 A charging adapter • Charging time may vary depending on charging adapter specifications, battery level, ambient temperature and charging conditions.

• Design and specifications are subject to change without notice.

Precautions for Use

Precautions before skin care

- It is not recommended to use the product over other parts of the body except face.
- When this product is used, there may be skin irritation such as blushing or tingling depending on the user's skin condition. If such irritation persists, stop using it immediately and consult a dermatologist before using it again.
- Those who have a facial skin disease or frequently suffer from facial dermatitis should consult a dermatologist before use.
- Those with photo-sensitivity such as sun allergy should not use the product.
- If insufficient skincare are used with the product, there may be skin irritation such as tingling.
- Using the product more frequently than the recommended use frequency may cause skin irritation.
(Intensive Care / Eye Care: 3 times a week, Cooling Care: Once a day)
- Using a cosmetic that contains the following ingredients may reduce its effectiveness, so be cautious of using the product.
(Example : Denatured alcohol, sunscreen, retinol, oil ingredients)
- Do not use the product while its head is fixed at one place, but keep moving it according to the guide in this user manual while using the product.
- Eye Care Head is recommended to be used around the eyes / mouth.
- Use Cooling Care only for skin care purposes.
- If a scratch or breakage occurs on the part of the head that comes into contact with the skin, stop using it and have it serviced at a service center.

Safety cautions

- **CAUTION concerning the Power Cord**
The Power Plug is the disconnecting device. In case of an emergency, the Power Plug must remain readily accessible.
- **CAUTION:** Only use the AC adapter supplied with this device. Do not use a power supply from another device or another manufacturer. Using any other power cable or power supply may cause damage to the device and void your warranty.
- When charging the product, be cautious not to get the charging part to be stained with foreign substances such as moisture or cosmetics. Otherwise, it may cause damage to the product or a fire.
 - **WARNING:** Do not use this appliance near bathtubs, showers, basins or other vessels containing water.



- Do not directly look into the LED light while the product is in operation.
- Cooling Care should not be done continuously. When Cooling Care is used continuously, it may shut down by itself after Voice Guide to maintain the proper temperature inside the product. In this case, wait for 60 minutes before using it again.

! For safety precautions, refer to the manual provided separately.


Setting Voice Guide



While the power of the Intensive / Eye Care is turned on, whenever both the Intensive / Eye Care MODE button and the Intensive / Eye Care LEVEL button are pressed for about 3 seconds simultaneously, the Voice Guide setting changes.

- Voice guide → Sound effects
→ Turn off the guide (Mute)

Setting Language

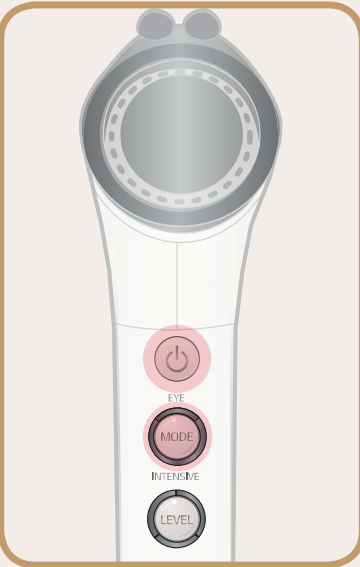


While the power of the Intensive / Eye Care is turned on, whenever both the Intensive / Eye Care Power button and the Intensive / Eye Care MODE button are pressed for about 3 seconds simultaneously, the used language changes.

Care Mode

Care Name	Head Used	Operational Status LED
Intensive Care (6 minutes)		
Eye Care (3 minutes)		
Cooling Care (3 minutes)		

Turning On/Off the Warming Mode in Intensive Care

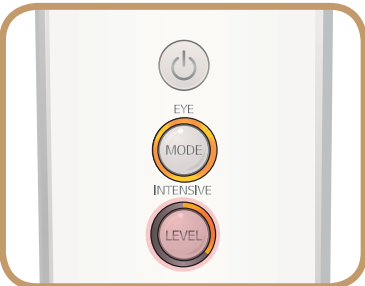


While the power of the Intensive / Eye Care is turned on, press the Intensive / Eye Care LEVEL button for about 3 seconds each and select the Warming Mode (on/off).

- The Warming Mode is not supported in the Eye Care and Cooling Care.

- ※ When the power of Intensive / Eye Care of the product is turned off and on, the Warming Mode is set to “On” automatically.
- ※ However, if the heat from the Warming Mode is too intense, use it while turning off the mode.
- ※ In the Warming Mode, you cannot change the temperature.


Intensive Care (For the Warming Mode, refer to Turning On/Off the Warming Mode in Intensive Care.)

Recommend to use 3 times a week It is recommended to use skincare that have been used previously or are suitable for your skin type. (Galvanic gel or vitamin serum, etc.)

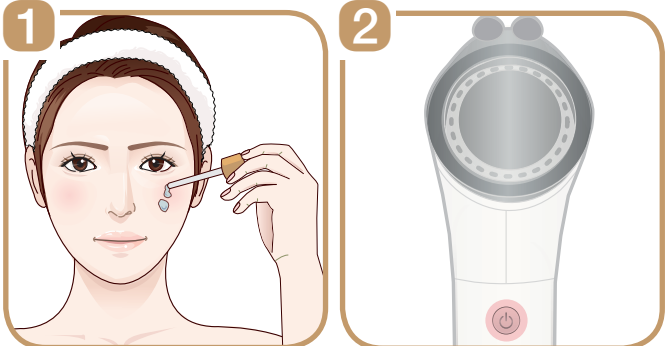

After washing your face, apply skincare sufficiently to each part to be cared for.

Press the Intensive / Eye Care Power button to turn on the product.

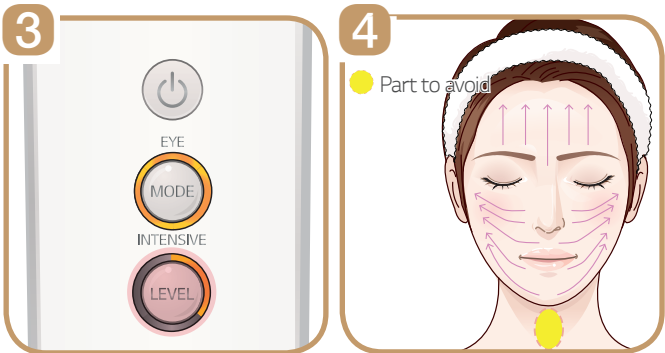


Press the Intensive / Eye Care LEVEL button to select the intensity of the Intensive Care (level 1 to 3) that suits you.

Massage the forehead / cheeks / chin areas for 6 minutes by lifting the Intensive Care Head slowly as shown in the picture.

- ※ Use Eye Care mode around your eyes.


Eye Care

Recommend to use 3 times a week It is recommended to use skincare that have been used previously or are suitable for your skin type. (Galvanic gel or vitamin serum, etc.)

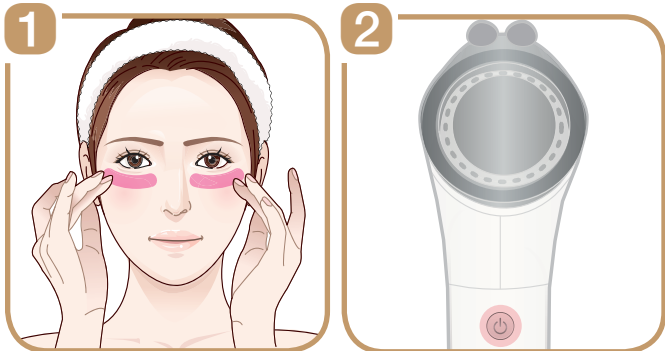

After washing your face, apply skincare sufficiently to the parts to be cared for.

Press the Intensive / Eye Care Power button to turn on the product.

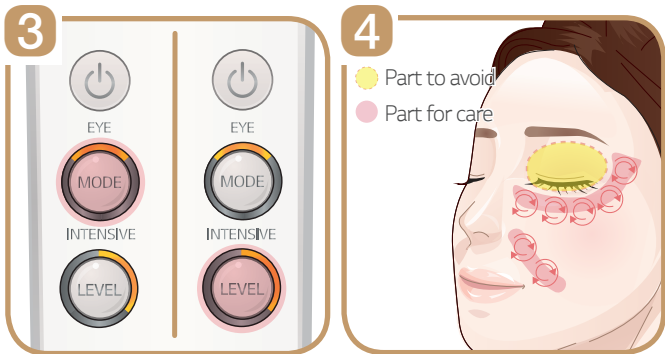


After selecting the Eye Care mode by pressing Intensive / Eye Care MODE button, select intensity (Level 1~3) by pressing Intensive / Eye Care LEVEL button.

As shown in the figure, massage the specific eye / mouth region with the eye care head slowly for 3 minutes.

Cooling Care

Recommended to use once a day, Cooling Care after Intensive Care / Eye Care helps to soothe the skin.

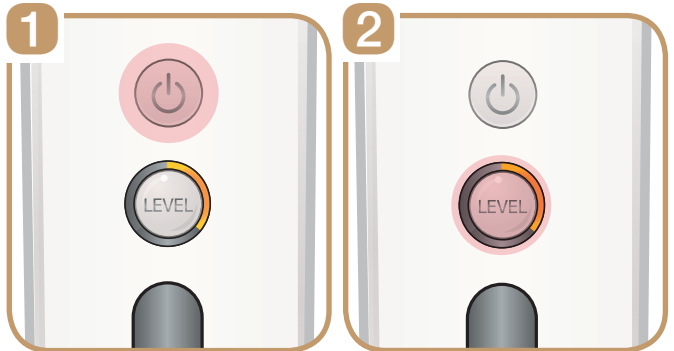

Press the Cooling Care Power button to turn on the product.


Press the Cooling Care LEVEL button to select the temperature of the Cooling Care (Level 1~3*) that suits you.

*The higher the level, the lower the temperature



- ※ Do not use Cooling Care continuously. When Cooling Care is used continuously, it may shut down by itself after Voice Guide to maintain the proper temperature inside the product. In this case, wait for 60 minutes before using it again.

Precaution for using Intensive Care mode

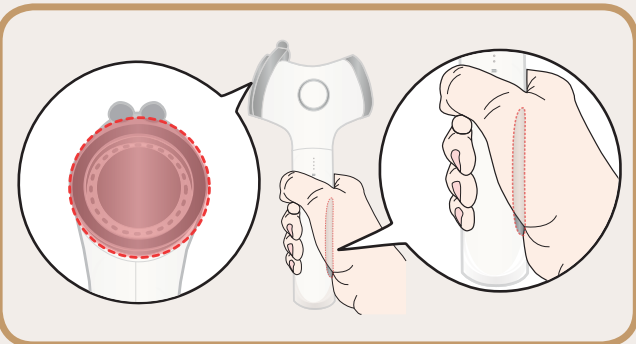

Gently wrap the skin contact electrode on the back of this product so that it touches the palm of your hand. At the same time, the product works only when the marked Intensive Care Head is in good contact with the skin of the face.

Q&A

Q How do I maintain the product?

After the care session, if foreign substances such as cosmetics remain on the product or on its head, wipe the part with a wet wipe or dry cloth, then dry the product in a cool and dry place.



Do not use strong chemicals such as alcohol, benzene, or acetone.



Be cautious for water or cosmetics not to enter the ventilation holes of the product.